

MADE IN FRANCE



# GENERATION OF **EXCELLENCE**

— FACE & BODY —



ADVANCE**EB**EAUTY

# THE SKIN'S SCIENCE RESULTS

Advance Beauty is a french company specialized in beauty treatment for face and body. Our scientific and technological know how provides innovative and performance solutions for complete treatment. In accordance with our standard quality organization, we fully apply our expertise to our customers' success.

The causes of skin ageing are known to stem from multiple and complex origins, that only a combination of specific and complementary treatments can actually set straight. Relying on more than 20 years of experience, Advance Beauty proposes a new strategy to erase the marks of time by rejuvenating the face and reshaping the body.

This strategy combines innovative treatments, stimulating the natural or physiological process thanks to our equipments and "external" treatment with the application of our ADM dermo-cosmetics, available for professional and home care users.

All our treatments are based on human physiology making them totally painless and comfortable for a perfect integration in any SPA center. Thanks to our DUOLIFT SEQUENTIAL system and to our easy and detailed **protocols**, you will benefit from an exclusive know how. We guarantee you visible, instant and long lasting results.



## THE RANGE DUOLIFT SEQUENTIAL FACE & BODY



### FRIENDLINESS

**A simple touch of the screen starts your selected treatment !**

- Extra large coloured HD screen 9" to guide your choice.
- Instant access to precise and exclusive protocols or to the Monitor program in order to customize your practice.
- Display of all of the treatment's parameters for maximum comfort.

### SAFETY

- Ultrasound probes equipped with a bio-compatible titanium resonator.
- RF temperature limiter: the sensor enables you to control and limit the temperature applied to the client's skin for comfort and absolute security.
- Compliance with CE norms.

### PROFITABILITY

The Sequential mode allows you to reduce the duration of the application and improve your profitability: a full facial in only 40 minutes compared to 1 hour 15 minutes using manual techniques.

The results are consistent and provide total satisfaction to your customers.



Advance Beauty goes always to the limits of technological feasibility to make the best use of known techniques because of their efficiency and safety.

1 | **SOFT PEELING**

A preliminary stage to start every protocol, the spatula **eliminates the dead cells, impurities and pollution**. It prepares the skin perfectly and provides its softness and radiance.

2 | **Computerized RADIOFREQUENCY**

Reshapes the deep layers of the epidermis to **skin redensify** providing long term effects.

3 | **SEQUENTIAL ULTRASOUND**

(Face & Body)  
Guarantee a **maximal reactivation of the metabolic processes to smooth tissues and skin breath**.  
The sonophoresis increases **the penetration** of dermo-cosmetic actives.

*So, facial cares will provide a face «lifting» effect thanks to skin rejuvenation and visible instant revitalisation.  
Ultrasound will also reactivate the metabolic lipolysis process in order to stimulate the loss in volume and reduce the fatty deposits («orange peel» appearance) during slimming treatments.*

4 | **Revolving LUMINOLOGY**

(with sequential ultrasound)  
Combined with our second generation of Luminology, the sequential system conveys in visible and **stimulating light** the evolution of the ultrasound waves throughout our protocols, indicating the constant changes to the stimulation process.  
This luminous energy **amplifies the targeted result**.



**FACE TREATMENT**



**BODY TREATMENT**



**AN EXPERTISE TO PROVIDE AND COMBINE THE BEST TECHNIQUES  
FOR SPECTACULAR RESULTS**



# HOW TO CHOOSE YOUR DUOLIFT SEQUENTIAL ?

Successors of our DUOLIFT Color, the DUOLIFT Sequential B can be easily carried out for home care.

## DUOLIFT SEQUENTIAL BV

Compact equipment (Box) for facial treatments, coming with 4 probes :

- 1 soft-peeling spatula,
- 2 sequential ultrasound face probes and Luminology,
- 1 radio frequency probe for face

## DUOLIFT SEQUENTIAL BCV

A standard with high performances for all your working needs: Face & Body equipped with avec 7 probes.

In addition to the facial features from the BV model, it comes with:

- 2 sequential ultrasound body probes and Luminology,
- 1 radio frequency probe for body

### TECHNICAL CHARACTERISTICS

Electricity supply	100-230 V AC / 50-60 Hz
Power	200 VA
Dimensions	33 (L) x 49(l) x 29 (h) cm
Weight	8 Kg



## DUOLIFT SEQUENTIAL RCV

A prestigious and complete unit for Face and Body. Coming with 7 probes, it has all functions from the BCV model in a luxurious presentation.

Mobile and design, it'll build up the image of your beauty center and fits beautifully into a SPA.



### TECHNICAL CHARACTERISTICS

Electricity supply	100-230 V AC / 50-60 Hz
Power	200 VA
Dimensions	33 (h) x 49(l) x 108 (h) cm
Weight	12 Kg



# PROBES SPECIFICATIONS

## ULTRA-VIBRATING SPATULA

Vibrating frequency	25 KHz nominal (piezo-electric activation)
Weight	210 g



## RF BODY PROBE

Frequency	450 KHz (nominal + overmodulation)
Power	40 W
Weight	210 g



## RF FACE PROBE

Frequency	450 KHz (nominal + overmodulation)
Power	40 W
Weight	120 g



## ULTRASOUND BODY PROBES

Frequency	1 MHz
Power maxi	1,5 W/cm <sup>2</sup>
Emission mode	Sequential
Active area	17 cm <sup>2</sup>
Sequential Luminology	included
Weight	450 g



## ULTRASOUND FACE PROBES

Frequency	2,8 MHz
Power maxi	1,5 W/cm <sup>2</sup>
Emission mode	Sequential
Titanium active area	0,8 cm <sup>2</sup>
Sequential Luminology	included
Weight	90 g



# SEQUENTIAL ULTRASOUND

Benefits from ultrasound waves have been known for decades. Advance Beauty was the pioneer in using this technique in facial and neck skincare.

Ultrasound generates an acoustic pressure to the skin and creates an internal micro massage which main biological effects are:

- The reactivation of micro-circulations (lymphatic and venous)
- The fibrolysis of tissue to smooth lines and wrinkles, or cellulite aspect.

Now, Advance Beauty innovates offering SEQUENTIAL ULTRASOUND that introduces new and exclusive benefits compared to existing practice:

- Automated treatment: the sequence of different stages is easy and intuitive to guarantee a simpler usage and to save time due to shorter sessions.

- Re-activate and maximize the metabolic process thanks to the cycle of sequential ultrasonic waves that generate continuous stimuli (internal micro massage), without any habituation and maintains the cells in a constant state of alert.

- The sonophoresis principle, that allows the penetration of active cosmetic ingredients through the cutaneous barrier, remains a major action of the sequential principle.

## SEQUENTIAL PRINCIPLE PATENTED SYSTEM

*In case of «alert» (physical reaction coming from external stimulations), the physiological processes will be reactivated, improving micro-circulations and the production of fibroblasts. In a classic use of ultrasound on the skin the «alert» process is very important during the first few minutes and through habituation loses its stimulating effects as the treatment progresses. The Sequential mode principle, patented by Advance Beauty, is based on a chain of ultrasonic sequences, different and non repetitive. So the stimulation remains constant throughout the treatment to increase the results and guarantee maximum efficiency for the 100% session's duration.*





# RADIOFREQUENCY

Radiofrequency (RF) is today the technique that cannot be ignored for significant non surgical skin tightening.

Radiofrequency heats the tissue of the skin internally via the energy coming from its high frequency current. Thanks to this localized temperature increase there are 2 principal actions on the fibroblasts:

- Collagen fibers contract, the deepest layers of the epidermis are reshaped.
- Increase of new collagen deposits restores the skins softness and tonicity over the weeks.

So we can get the intended firming effect in facial treatments as well as smoother skin usefull for body treatments.


On top of that, the localized heating produces hyperhemia. As a consequence, a strong re-hydration comes instantly to the superficial skin layers.

This water brought up to the superficial layers takes part in the «fleshy» skin effect and will also amplify the ultrasound results.



Fibroblasts are cells that secrete proteins creating connective tissue fibers. These cells maintain skin structure by constantly producing new Collagen and other structural proteins. Effective heating of skin tissue does two things in regards to fibroblasts:

- Existing collagen fibers contract,
- Fibroblasts' metabolism is restarted and accelerated, creating new collagen deposits.

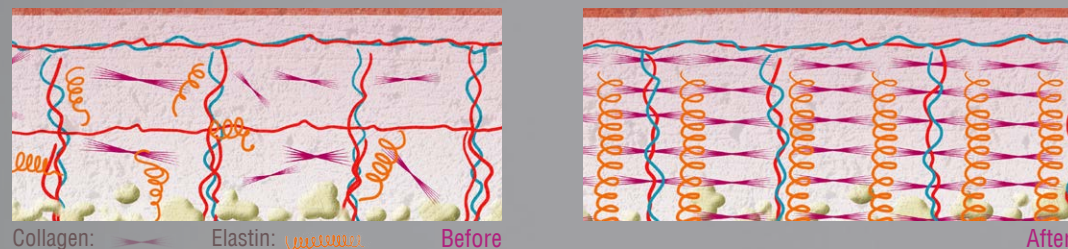

# ADM SKIN SCIENCE 3.0

NEW ANTI-AGING AND REPAIRING LINE

Highly concentrated, the complexes associated in each of our formulas help us to fight the causes and consequences of skin ageing and marks of time.

These formulas have been carefully selected to increase the performance of our protocols in synergy with the DUOLIFT SEQUENTIAL.

A short and coherent range of 6 references designed for professional use and/or retail.



**RADIANCE PEELING LOTION**  
**LIFT SERUM**  
**RESTRUCTURING**  
**ACNO-LYSE**  
**GEL COLLAGEN**  
**DERMA-LYSE**



ADVANCE **E** BEAUTY



**ADVANCE**  
**BEAUTY**

Rovaltin, parc 45 parallèle  
5 Rue Olivier De Serres  
26958 Valence Cedex 9  
FRANCE

Tél : +33 (0)4 75 48 17 54  
Fax : +33 (0)4 75 47 74 93

[Info@advance-beauty.com](mailto:Info@advance-beauty.com)  
[www.advance-beauty.com](http://www.advance-beauty.com)

