



MADE IN FRANCE

ADVANCE **EB**EAUTY

NEWS SEQUENTIAL

SLIMMING STATION
SHOOT RIGHT YOUR TARGET !



Based on a 20 years old expertise and in line with its combination of techniques philosophy Advance Beauty invents the NEWS SEQUENTIAL.

The only equipment that enable to simultaneously and selectively treat superficial, deep and localized fat thanks to exclusive and non-invasive technologies.

The NEWS SEQUENTIAL is easy to use thanks to many details that offer an ergonomic application. The independent operator treatments associated to the localized actions of the two hand probes ensure efficiency with a high profitability.

The NEWS SEQUENTIAL is the inescapable tool for getting success in all slimming programs.



SEQUENTIAL ULTRASOUND

Innovative technology patented by Advance Beauty

The benefits of ultrasound waves are known for decades.

Advance Beauty was the pioneer in using this technique in the anti-ageing and slimming treatments.

The ultrasound waves are super-fast mechanical vibrations which generate an acoustic pressure to the skin. They create an internal micro massage which main biologicals effects are:

- The reactivation of micro-circulations (lymphatic and venous)
- The fibrolysis of tissue to smooth the cellulite aspect.

Now, Advance Beauty innovates offering SEQUENTIAL ULTRASOUND that introduces new and exclusive benefits:

- Automated treatment: the sequence of different stages is easy and intuitive to guarantee a simpler usage and to save time due to shorter sessions.
- Re-activate and maximize the vital processes of the skin: the cycle of sequential ultrasonic waves generates continuous stimuli without any habituation and maintains the cells in a constant state of alert.
- The sonophoresis principle, that allows the penetration of active cosmetic ingredients through the cutaneous barrier, remains a major action of the sequential principle.

THE SEQUENTIAL PRINCIPLE

Patented system

In case of « alert » (physical reaction coming from external stimulations), the physiological processes will be reactivated, improving the lipolysis and the micro-circulations.

In a classic use of ultrasound on the skin the « alert » process is very important during the first few minutes and through habituation loses its stimulating effects as the treatment progresses.

The Sequential mode principle, patented by Advance Beauty, is based on a chain of ultrasonic sequences, different and non-repetitive. So the state of “alert” remains constant throughout the treatment to increase the results and guarantee maximum efficiency for the full duration of the session.



4 COMPLEMENTARY ACTIONS

1 | IMMEDIATE REDUCTION OF LOCALIZED FAT DEPOSITS – VISIBLE RESULTS

The LOCA-LYSE probe is manually used on limited zones, having strong excess weight, in order to reduce remarkably the adipose tissue.



2 | LIPOLYSIS OF SUPERFICIAL FAT

The 10 ultrasound pads of low frequency act in a mechanical way on superficial fat by activating the process of metabolic lipolysis.

Thanks to their 30 impact points, they participate in the reduction of the orange peel appearance and authorize the treatment of a very vast surface in a single session.



3 | LIPOLYSIS OF DEEP FAT


The 2 ultrasound pads of middle frequency stimulate the nerve endings of sympathetic nervous system, which by induction, favor the lipolysis of the lower layers and thus the destocking of deep fat.



4 | FIRMING AND SKIN SMOOTHING

The radio frequency probe allows to obtain a significant strengthening without surgery.

It has become today a key technique.

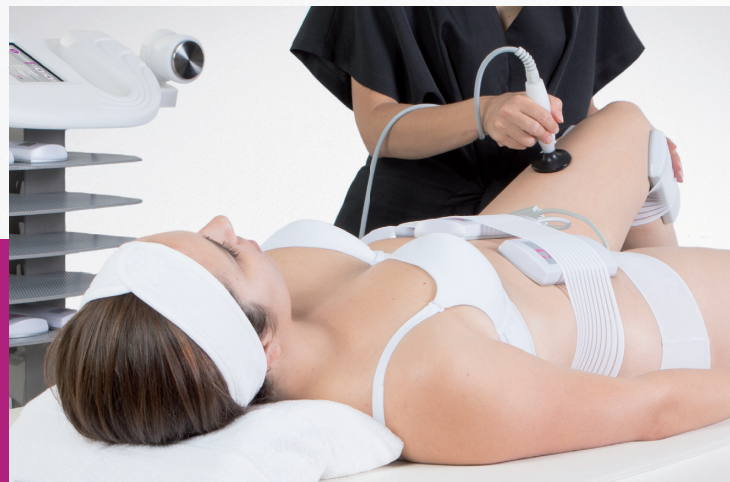

Thus, the NEWS SEQUENTIAL realize a complete slimming action associated with a reshaping effect.

PROVEN RESULTS

Clinical study made by Pr. Humbert – CHU Besançon
France in 2016 on 20 women.

Study with independent operator plates only.

**On 90 % of the study subjects we obtained the tightening effect
a reduction of the orange skin peel appearance.**

A significant clinical improvement of the cellulite has been found.

Before



After




Study realized on 28 volunteers
with the LOCA-LYSE probe only.

AFTER 3 SESSIONS
A significant circumference loss

- 3.88 CM
ON THE ABDOMEN

- 4.35 CM
ON THE HIPS

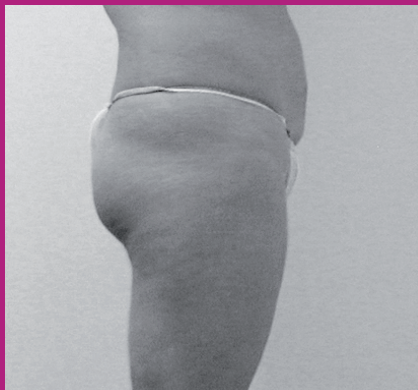

TRANSDUCERS HIGH PERFORMANCE

The NEWS SEQUENTIAL generates premium quality sequential ultrasound thanks to exclusive ceramics which provide a very high efficiency and a real balanced distribution to the whole area of the transducers (probes and pads). This EXCLUSIVE advanced technology guarantees the best performance.

1

5 mobile shelves to apply the gel and organize very easily your session: saving of time and space.

2


The RF probe comes with an integrated control and regulation system of the temperature insuring comfort and safety of the client.

3

The LOCA-LYSE probe has a resonator in biocompatible titanium which ensures perfectly the waves transmission.

2

3

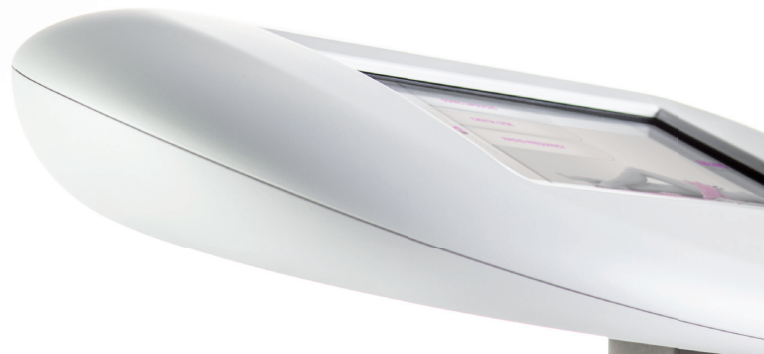

DESIGN, USER-FRIENDLINESS AND ERGONOMY

A simple contact on the screen and the chosen treatment begins!

- Extra-large color screen HD to guide your selection in an intuitive way.
- Display of all the parameters of the treatment for an optimal comfort.

The NEWS Sequential is perfectly movable and its congestion optimized.
The storage system of the various applicators offers you a practical access and only 5 cords connect the device to the 12 pads.



TECHNICAL SPECIFICATIONS

Power supply	100-230 V AC / 50-60 Hz
Emission type	sequential, continuous or modulated
Maximum ultrasound power	2 W/cm²
Display	8" colored
Dimensions (cm)	110 (h) x 51(l) x 30(p)
Weight without resonators	30 Kg
Storage plates system	5 swivel shelves included
Connection to the system	USB 2 type

ACCESSORIES CHARACTERISTICS

Ultrasound L.F pads	Nominal frequency	30 KHz
	Number of cells / pad	3
	Number of pack of 3 pads	2
	Number of pack of 2 pads	2
	Total number of cells	30
	Active area / plate	110 cm² (17 x 6,5 cm)
	Operation indicator	integrated leds
	Unit weight	230 g
Ultrasound M.F pads	Nominal frequency	500 KHz
	Number of cells / pad	3
	Number of pack of 2 pads	1
	Total number of cells	6
	Active area / plate	110 cm² (17 x 6,5 cm)
	Operation indicator	integrated leds
	Unit weight	390 g
LOCA-LYSE probe	Frequency	40 KHz
	Maximum power	25 W
	Active area	24 cm²
	Unit weight	700 g
RF probe	Type	Resistive
	Frequency	450 KHz
	Maximum power	70 W
	Active area	20 cm²
	Unit weight	210 g

2 sets of elastic belts of various lengths are delivered with the device to fix the plates on the skin.



The NEWS SEQUENTIAL is unique to propose such a synergy of techniques in slimming and cellulite treatments. NON-INVASIVE treatment.

ADVANCE  BEAUTY

Rovaltain, parc 45 parallèle
5 Rue Olivier De Serres
26958 Valence Cedex 9
FRANCE

Tél : +33 (0)4 75 48 17 54
Fax : +33 (0)4 75 47 74 93

Info@advance-beauty.com
www.advance-beauty.com

